



NITINOL GUIDEWIRES

Threefold benefit through layer technology



THREEFOLD BENEFIT THROUGH LAYER TECHNOLOGY


Under the RÜSCH brand name, Teleflex provides well-known quality products for almost every area of urology – from transurethral drainage to nephrostomy. More than 100 years of know-how have gone into each of our products, and we have been able to exploit this probably unique experience to make them as simple, comfortable and safe as possible, for patients, doctors and other medical staff.

However, we are not content with what we have achieved. Our constant search for perfection through continuous research and development plus the use of new technologies has resulted in medical products, the quality of which sets new standards throughout the world.

Our comprehensive range of hydrophilic coated high-quality nitinol guidewires is suitable for a wide range of urological applications such as the inspection of hollow tubular organs or the insertion of medical devices.


The guidewires can be placed in the efferent urinary tract, or used in connection with suprapubic and percutaneous puncture channels.



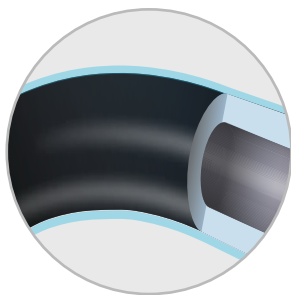


BLACK POLYURETHANE JACKET

makes the guidewire easily visible in dark endoscopy rooms. The highly radiopaque material enables control of the position using x-ray images.

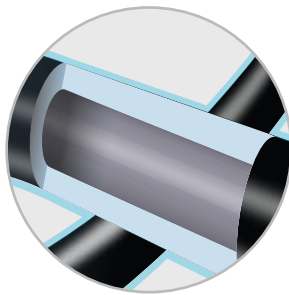

NITINOL CORE

makes the guidewire extremely resistant to kinking and gives it high rotational stability, consequently enabling precise handling during treatment.

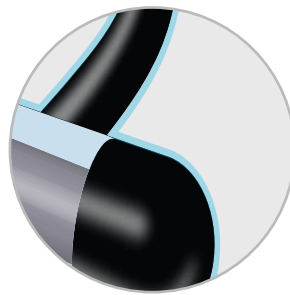


HYDROPHILIC COATING

provides fast and homogeneous hydration and significantly improves sliding properties. This facilitates the application, especially when faced with narrow and obstructed passages.

KEY BENEFITS

- excellent resistance to kinking and high rotational stability
- high visibility in dark endoscopy rooms
- high radiopacity
- excellent sliding properties
- atraumatic application
- broad range of guidewires
- comes in a dispenser with Luer-lock adapter and introducer
- sterile, single use
- latex-free

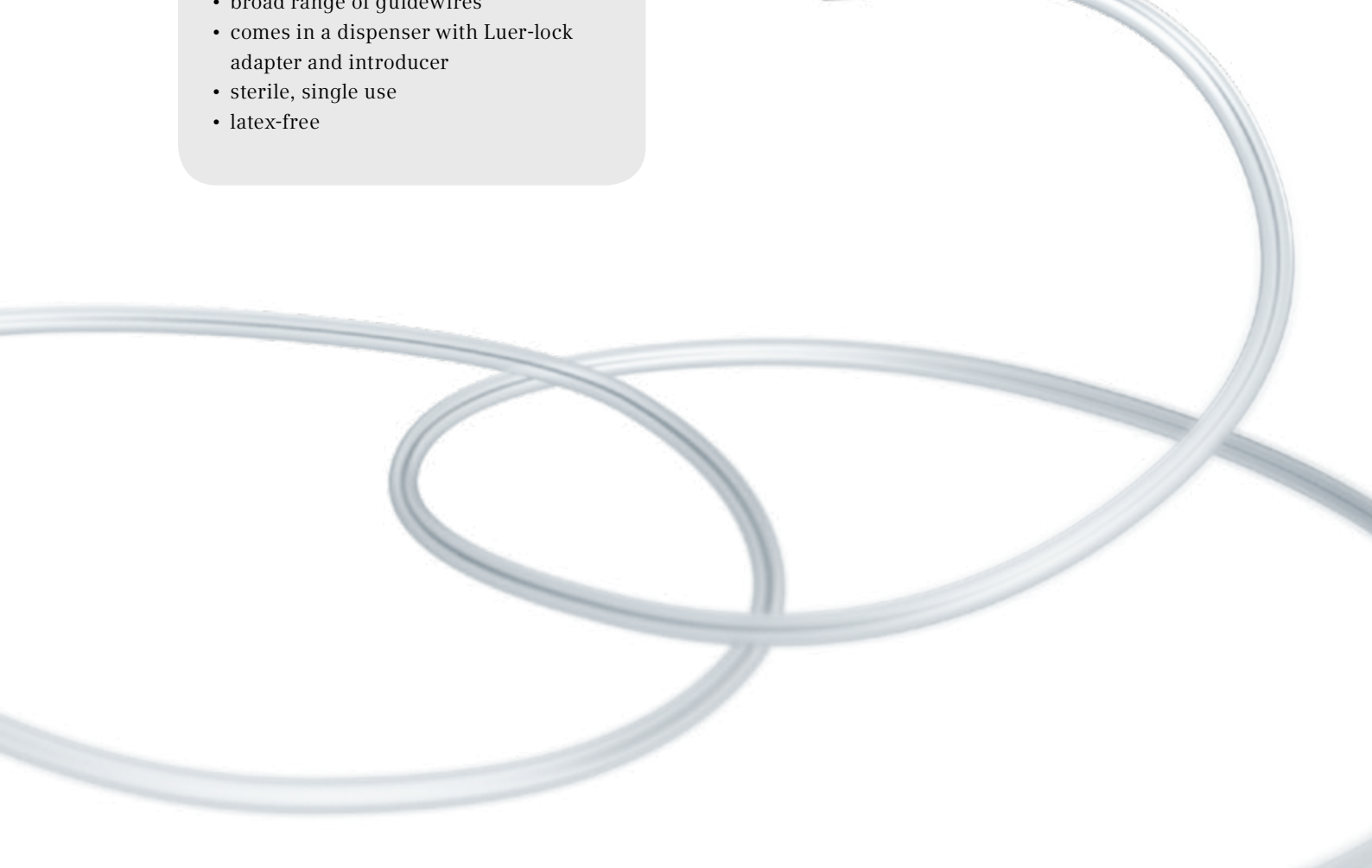

SUPERIOR SAFETY

All our nitinol guidewires, whether with a flexible or a stiff shaft, are provided with one flexible safety tip which protects the sensitive mucosa of the urinary tract from irritation and lesions.

The transition zone between the shaft and the flexible safety tip is particularly smooth to preserve the delicate urinary tract or the puncture channel.

THE CORRECT GUIDEWIRE FOR EVERY APPLICATION

Depending on the type of application and the devices the guidewire is to be used with, you can choose between those with a flexible or a stiff shaft, each available in different lengths and different sizes. The guidewires are available either with a straight or an angled flexible safety tip.

For special applications, we can also provide flexible and stiff nitinol guidewires with a longer and extra soft angled flexible safety tip.

COMFORTABLE APPLICATION


All our nitinol guidewires come in a dispenser with Luer-lock adapter for easy activation of the hydrophilic coating by means of a Luer or Luer-lock syringe.

Additionally, the dispenser features an introducer to facilitate insertion of the guidewire's flexible safety tip.





ORDER INFORMATION

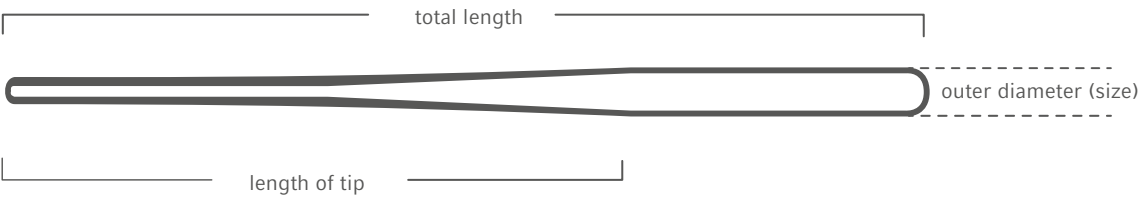

FLEXIBLE NITINOL GUIDEWIRES
WITH STRAIGHT STANDARD TIP (13 CM)

	REF.	LENGTH	SIZE	QTY
	680160-150018	150 cm	0.018"	5
	680160-150035	150 cm	0.035"	5
	680160-150038	150 cm	0.038"	5
	680160-180018	180 cm	0.018"	5
	680160-180035	180 cm	0.035"	5
	680160-180038	180 cm	0.038"	5
	680160-260035	260 cm	0.035"	5
	680160-260038	260 cm	0.038"	5

FLEXIBLE NITINOL GUIDEWIRES
WITH ANGLED STANDARD TIP (13 CM)

	REF.	LENGTH	SIZE	QTY
	680161-150018	150 cm	0.018"	5
	680161-150035	150 cm	0.035"	5
	680161-150038	150 cm	0.038"	5
	680161-180018	180 cm	0.018"	5
	680161-180035	180 cm	0.035"	5
	680161-180038	180 cm	0.038"	5
	680161-260035	260 cm	0.035"	5
	680161-260038	260 cm	0.038"	5

FLEXIBLE NITINOL GUIDEWIRES
WITH SOFT ANGLED TIP (19 CM)

	REF.	LENGTH	SIZE	QTY
	680162-150035	150 cm	0.035"	5
	680162-150038	150 cm	0.038"	5

STIFF NITINOL GUIDEWIRES
WITH ANGLED STANDARD TIP (13 CM)

	REF.	LENGTH	SIZE	QTY
	680170-150035	150 cm	0.035"	5
	680170-150038	150 cm	0.038"	5
	680170-180035	180 cm	0.035"	5
	680170-180038	180 cm	0.038"	5

STIFF NITINOL GUIDEWIRES
WITH ANGLED STANDARD TIP (13 CM)

	REF.	LENGTH	SIZE	QTY
	680171-150035	150 cm	0.035"	5
	680171-150038	150 cm	0.038"	5
	680171-180035	180 cm	0.035"	5
	680171-180038	180 cm	0.038"	5
	680171-260035	260 cm	0.035"	5
	680171-260038	260 cm	0.038"	5

STIFF NITINOL GUIDEWIRE
WITH SOFT ANGLED TIP (19 CM)

	REF.	LENGTH	SIZE	QTY
	680172-150035	150 cm	0.035"	5

Teleflex is a leading global provider of specialty medical devices used for diagnostic and therapeutic procedures in critical care, urology and surgery. Our mission is to provide solutions that enable healthcare providers to improve outcomes and enhance patient and provider safety. We specialise in devices for general and regional anaesthesia, cardiac care, respiratory care, urology, vascular access and surgery and we serve healthcare providers in more than 150 countries. Teleflex also provides specialty products or medical device manufacturers.

Our well known brands include ARROW®, DEKNATEL®, GIBECK®, HUDSON RCI®, KMEDIC®, LMA™, PILLING®, PLEUR-EVAC®, RÜSCH®, SHERIDAN®, TAUT®, TFX OEM®, VASONOVA™, VIDACARE® and WECK®, all of which are trademarks or registered trademarks of Teleflex Incorporated.

Teleflex global operations: Australia, Austria, Belgium, Canada, China, Czech Republic, France, Germany, Greece, India, Ireland, Italy, Japan, Malaysia, Mexico, Netherlands, New Zealand, Portugal, Singapore, Slovak Republic, South Africa, Spain, Switzerland, United Kingdom, Uruguay and USA.

YOUR INTERNATIONAL CONTACTS:

TELEFLEX HEADQUARTERS INTERNATIONAL, IRELAND

Teleflex Medical Europe Ltd., IDA Business and Technology Park,
Dublin Road, Athlone, Co Westmeath
Phone +353 (0)9 06 46 08 00 · Fax +353 (0)14 37 07 73
orders.intl@teleflex.com

AUSTRALIA/NEW ZEALAND +61 (0)3 9081 0600

AUSTRIA +43 (0)1 402 47 72

BELGIUM +32 (0)2 333 24 60

CHINA (SHANGHAI) +86 (0)21 6163 0965

CHINA (BEIJING) +86 (0)10 6418 5699

CZECH REPUBLIC +420 (0)495 759 111

FRANCE +33 (0)5 62 18 79 40

GERMANY +49 (0)7151 406 0

GREECE +30 210 67 77 717

INDIA +91 (0)44-2836 5040

ITALY +39 0362 58 911

JAPAN +81 (0)3 3379 1511

NETHERLANDS +31 (0)88 00 215 00

PORTUGAL +351 22 541 90 85

SINGAPORE +65 6439 3000

SLOVAK REPUBLIC +421 (0)3377 254 28

SOUTH AFRICA +27 (0)11 807 4887

SPAIN +34 918 300 451

SWITZERLAND +41 (0)31 818 40 90

UNITED KINGDOM +44 (0)1494 53 27 61

For detailed information see www.teleflex.com

The products in this catalogue may not be available in all countries. Please contact your local representative. All data current at time of printing (08/2014). Subject to technical changes without further notice.

94 03 49 - 00 00 01 · REV A · MC / WM · 08 14 01