

## RÜSCH U-BAG

LATEX-FREE CLOSED URINARY DRAINAGE SYSTEM WITH CHAMBER FOR HOURLY MEASUREMENT



### SAFE AND HYGIENIC: MAXIMUM SAFETY FOR THE PATIENT AGAINST NOSOCOMIAL URINARY TRACT INFECTIONS

About 90% of all urinary tract infections are caused by urinary drainage. Apart from the health risks for the patient himself, such infections also present a big economic problem. Rüsche offers a variety of products which reduce the risk of infection.

The closed urinary drainage system represents one of the most important elements of prophylaxis against infections. Rüsche urinary drainage systems meet all the hygienic demands which are made on a high-quality closed system for long-term medical care.


### Optimum pressure balance and excellent flow performance

Double bacteria-proof and liquid-tight ventilation: at the drip chamber and the urine collector. Prevents the congestion of urine in the tube and the associated risk of infection.



### AT A GLANCE

- Urine port without needle
- Drip chamber with integrated check valve
- 400 ml mounted pivoting double measurement chamber with integrated 50 ml measurement chamber for small amounts
- Double bacteria-proof and liquid-tight ventilation
- Clear graduation in ml
- Flexible and non-kinking tube
- Optimum fastening of the system
- Non-dripping bottom outlet
- Latex-free

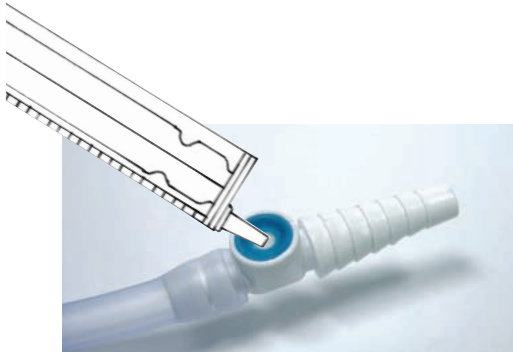

### Urine sampling without needle

Urine sampling is carried out at the ridged universal connector with a Luer syringe without a needle. Thus, needle stick injuries are eliminated and the risk of invading bacteria is greatly reduced. The connector is tapered opposite to the direction of extraction, thus preventing it from slipping out of the catheter funnel.



**Mounted 400 ml double measurement chamber, which can be pivoted for drainage** with integrated 50 ml measurement chamber for small amounts (5 ml graduation), with additional sampling port for fresh urine (1 ml graduation)



### Hygienic and quick emptying

Non-dripping bottom outlet for one-hand operation. Doghouse holder to fix the outlet tube, to prevent contamination, and to achieve additional floor clearance.



**Flexible tube, resistant to kinking** with sliding clamp and sheet clip

**2000 ml urine collector** vented, with clear graduation



**Vented drip chamber** with integrated check valve



**Optimum fastening** with velcro band and D-ring

## PRODUCT INFORMATION

Latex-free closed urinary drainage system with chamber for hourly measurement

Rüsch U-Bag			
Ref.	Capacity	Length of tube	Qty
67 30 64	2000 ml	120 cm	10
Sterile, for single use, latex-free			

Distributed by:

**GERMANY** • Teleflex Medical GmbH • P.O. Box 11 80 • 71385 Kernen  
 Phone +49 (0)7151 406 0 • Fax +49 (0)7151 406 150 • info.de@teleflexmedical.com  
**UNITED KINGDOM** • Teleflex Medical • High Wycombe • Bucks. HP12 3ST  
 Phone +44 (0)1494 53 27 61 • Fax +44 (0)1494 52 46 50 • info.uk@teleflexmedical.com

www.ruesch.de • www.teleflex.com