

Teleflex is a leading global provider of specialty medical devices used for diagnostic and therapeutic procedures in critical care, urology and surgery. Our mission is to provide solutions that enable healthcare providers to improve outcomes and enhance patient and provider safety.

We specialise in devices for general and regional anaesthesia, cardiac care, respiratory care, urology, vascular access and surgery and we serve healthcare providers in more than 130 countries. Teleflex also provides specialty products for medical device manufacturers.

Our well known brands include ARROW®, BEERE MEDICAL®, DEKNATEL®, GIBECK®, HUDSON RCI®, KMEDIC®, PILLING®, PLEUR-EVAC®, RUSCH®, SHERIDAN®, SMD®, TAUT®, TFX OEM®, VASONOVA™ and WECK®, all of which are trademarks or registered trademarks of Teleflex Incorporated.

Teleflex global operations: Austria, Belgium, Canada, China, Czech Republic, France, Germany, Greece, India, Ireland, Italy, Japan, Mexico, Netherlands, Portugal, Singapore, Slovak Republic, South Africa, Spain, Switzerland, United Kingdom, Uruguay and USA.

YOUR CONTACTS FOR EUROPE, THE MIDDLE EAST AND AFRICA (EMEA):

TELEFLEX MEDICAL HEADQUARTERS EMEA, IRELAND

Teleflex Medical Europe Ltd., IDA Business Park, Athlone, Co Westmeath
Phone +353 (0)9 06 46 08 00 · Fax +353 (0)14 37 07 73
orders.intl@teleflex.com

AUSTRIA +43 (0)1 402 47 72

BELGIUM +32 (0)2 333 24 60

CZECH REPUBLIC +420 (0)495 759 111

FRANCE +33 (0)5 62 18 79 40

GERMANY +49 (0)7151 406 0

GREECE +30 210 67 77 717

ITALY +39 0362 58 911

NETHERLANDS +31 (0)88 00 215 00

PORTUGAL +351 22 541 90 85

SLOVAK REPUBLIC +421 (0)3377 254 28

SOUTH AFRICA +27 (0)11 807 4887

SPAIN +34 918 300 451

SWITZERLAND +41 (0)31 818 40 90

UNITED KINGDOM +44 (0)1494 53 27 61

For detailed information see www.teleflex.com

The products in this catalogue are only available in the EMEA (Europe, Middle East, Africa) region. For further information contact your local representative. All data current at time of printing (06/2011). Subject to technical changes without further notice.

94 06 24 - 00 00 01 · REV A · MC / WM · 07 11 01



RÜSCH URINE BAGS


Simple, safe & hygienic



RIDGED UNIVERSAL CONNECTOR WITH NEEDLE-FREE URINE SAMPLING PORT

The sampling port helps to avoid needle-stick injuries and reduces the risk of infection. For taking a urine sample, it is advisable, after disinfection of the urine sampling port, to press the tip of a Luer slip syringe without needle into the centre of the rubber stopper of the sampling port, turning it slightly, until it penetrates.


Furthermore all our urine bags with needle-free sampling port are equipped with a specially designed ridged universal connector with protective cap – suitable for all common catheters. The connector is tapered opposite to the direction of catheter extraction, thus preventing it from slipping out of the catheter funnel.

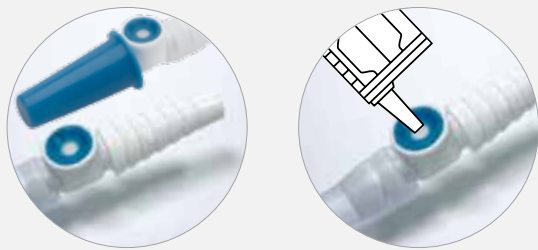

NON-DRIP OUTLET TAP

The outlet taps of our bags are designed to be user friendly. They are suitable for single handed use and for individuals with compromised dexterity.

Our leg bags are supplied with an integral extension tubing at the end of the outlet tap to facilitate additional connections.

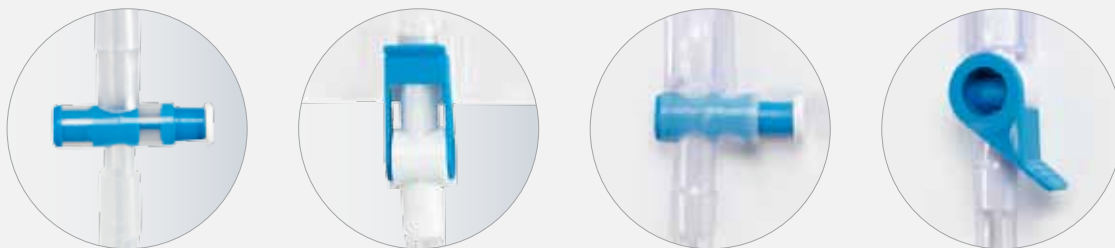

From left to right: slide tap (leg bag), lever tap (leg bag), slide tap (night bag) & lever tap (night bag)

LEG URINE BAGS

500 ML LEG URINE BAG WITH SLIDE TAP

- ridged universal connector with protective cap and needle-free urine sampling port
- flexible, kink-resistant inlet tubing
- non-reflux valve
- clear graduation in ml
- slide tap with prolongation tubing for further urine bags
- double packaging (each bag with instructions for use)
- 5 pairs of Velcro straps per packaging unit
- latex-free
- sterile, single use

500 ML LEG URINE BAG WITH SLIDE TAP				RUSCH
	REF.	CAPACITY	TUBE LENGTH	QTY
	673075-000500	500ml	10cm	10
	673076-000500	500ml	30cm	10

500 ML LEG URINE BAG WITH LEVER TAP

- ridged universal connector with protective cap and needle-free urine sampling port
- flexible, kink-resistant inlet tubing
- non-reflux valve
- clear graduation in ml
- lever tap with prolongation tubing for further urine bags
- double packaging (each bag with instructions for use)
- 5 pairs of Velcro straps per packaging unit
- latex-free
- sterile, single use

500 ML LEG URINE BAG WITH LEVER TAP				RUSCH
	REF.	CAPACITY	TUBE LENGTH	QTY
	673077-000500	500ml	10cm	10
	673078-000500	500ml	30cm	10

VELCRO STRAPS

Each packaging unit comes with 5 pairs of Velcro straps with silicone coating on the inside to ensure the leg bag is securely fixed.

The Velcro straps can also be ordered separately.

VELCRO STRAPS FOR LEG BAGS		
	REF.	QTY
	850360-420025	2

NIGHT URINE BAGS

2,000 ML NIGHT BAG WITH SLIDE TAP

- ridged universal connector with protective cap and needle-free urine sampling port
- flexible, kink-resistant inlet tubing
- robust double hook for secure fastening
- small drip chamber with integrated non-reflux valve
- bacteria-proof and liquid-tight ventilation at the bag
- clear graduation in ml
- slide tap with outlet housing for protection from contamination
- double packaging (each bag with instructions for use)
- latex-free
- sterile, single use

2,000 ML NIGHT BAG WITH SLIDE TAP				RUSCH
REF.	CAPACITY	TUBE LENGTH	QTY	
673072-000000	2,000 ml	120 cm	10	



2,000 ML NIGHT BAG WITH LEVER TAP

- ridged universal connector with protective cap and needle-free urine sampling port
- flexible, kink-resistant inlet tubing
- robust double hook for secure fastening
- small drip chamber with integrated non-reflux valve
- bacteria-proof and liquid-tight ventilation at the bag
- clear graduation in ml
- lever tap with outlet housing for protection from contamination
- double packaging (each bag with instructions for use)
- latex-free
- sterile, single use

2,000 ML NIGHT BAG WITH LEVER TAP				RUSCH
REF.	CAPACITY	TUBE LENGTH	QTY	
673073-000000	2,000 ml	120 cm	10	



BACTERIA-PROOF AND LIQUID-TIGHT VENTILATION AT THE BAG

ensures optimum pressure balance in the bag and reliable flow performance, even under high pressure load. Congestions of urine inside the inlet tubing are avoided and the associated risk of infection is reduced.



OUTLET HOUSING TO FIX THE OUTLET TUBING

ensures secure fixation of the outlet tubing, prevents potential contamination of the surrounding area and achieves additional floor clearance.



TELEFLEX RÜSCH URINE BAGS: SIMPLE, SAFE AND HYGIENIC

Trusted brands make Teleflex a reliable and strong partner. Built on a solid tradition of innovation, Teleflex is a global leader in superior medical supplies designed to help providers minimise risk and maximise outcomes for their patients. Our understanding of the importance to our customers of a full range of products has led to the development of a unique line of medical devices, all of which complement one another.

It is widely accepted that a large percentage of urinary tract infections are caused by contaminated drainage systems. For this reason, we offer a choice of urine collection systems in order that this risk may be minimised.

We offer high quality and easy to use urine bags for day and night. Our urine bags meet all the hygienic demands whilst remaining economical with no compromise to quality.

The RÜSCH urine drainage bags shown in this brochure have three features in common: they all have needle-free sampling ports, are latex-free and have a non-drip outlet tap.

You can find further details and the technical specifications of our products within this brochure.

**TELEFLEX – HIGH QUALITY MEDICAL
SUPPLIES FROM A SINGLE SOURCE**

